



VitalJacket® SDK


VitalJacket SDK v1.0.07 – Android BioLib library

LEGAL NOTICE AND DISCLAIMER

ATTENTION: Although **VitalJacket** is a certified medical device, its developer version is **NOT** certified for diagnosis usage. It is intended for R&D and development purposes **only**. Users of VJ SDK can submit their final developments to medical certification. All contents of our product are compliant with the European Medical Device directive 93/42/EEC but, being a developer's version, it's not certified.

Index

1. Import library:	4
2. Declare the Bluetooth permission(s) in your application manifest file:	4
3. In application properties select Android project build target:	4
4. Add BioLibSDK.jar to your project:.....	4
5. In manifest file, insert:	5
6. BioLib.SDK.jar	5
7. Print screens of Android application:.....	10
8. References.....	11
9. Control versions	12

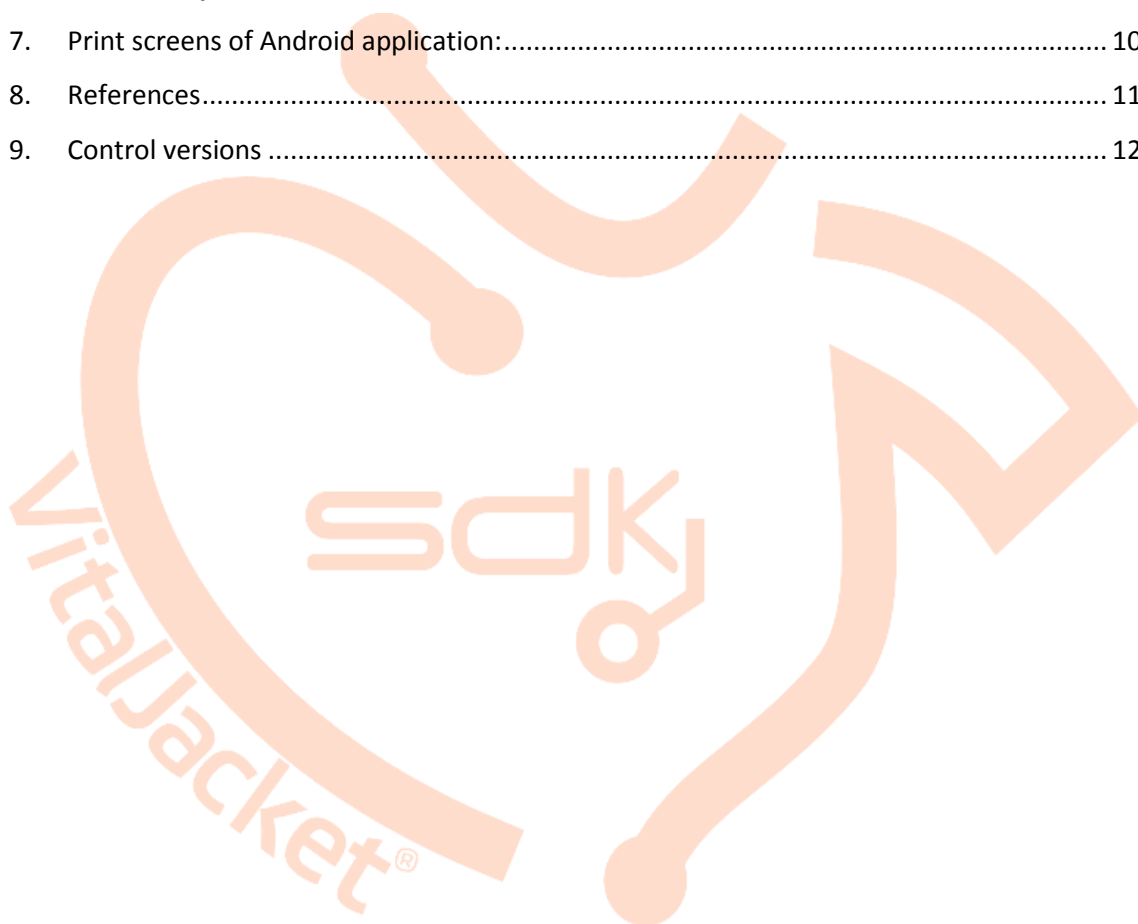


1. Import library:

```
import Bio.Library.namespace.BioLib;
```


2. Declare the Bluetooth permission(s) in your application manifest file:

```
<uses-permission android:name="android.permission.BLUETOOTH_ADMIN" />
<uses-permission android:name="android.permission.BLUETOOTH" />
```

3. In application properties select Android project build target:

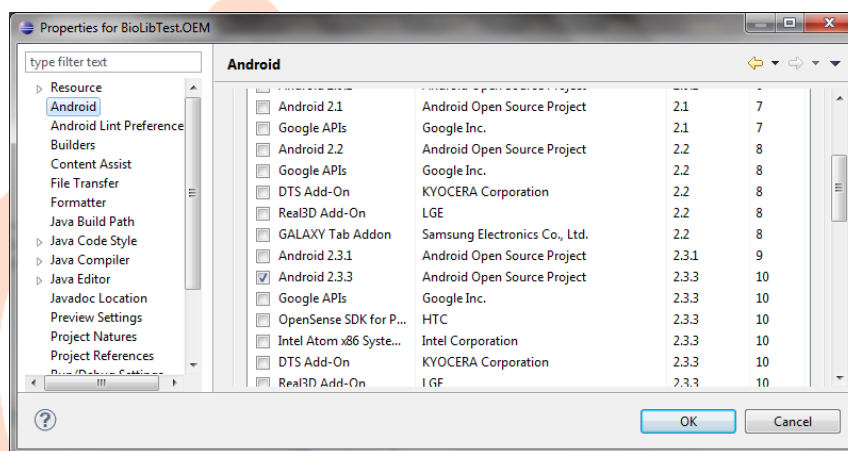



4. Add BioLibSDK.jar to your project:

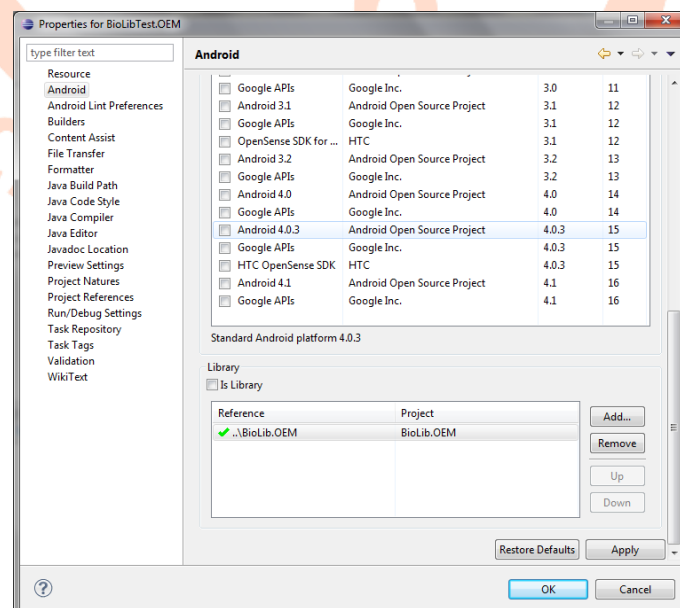

5. In manifest file, insert:

```
<uses-sdk android:minSdkVersion="10" />
```

6. BioLib.SDK.jar

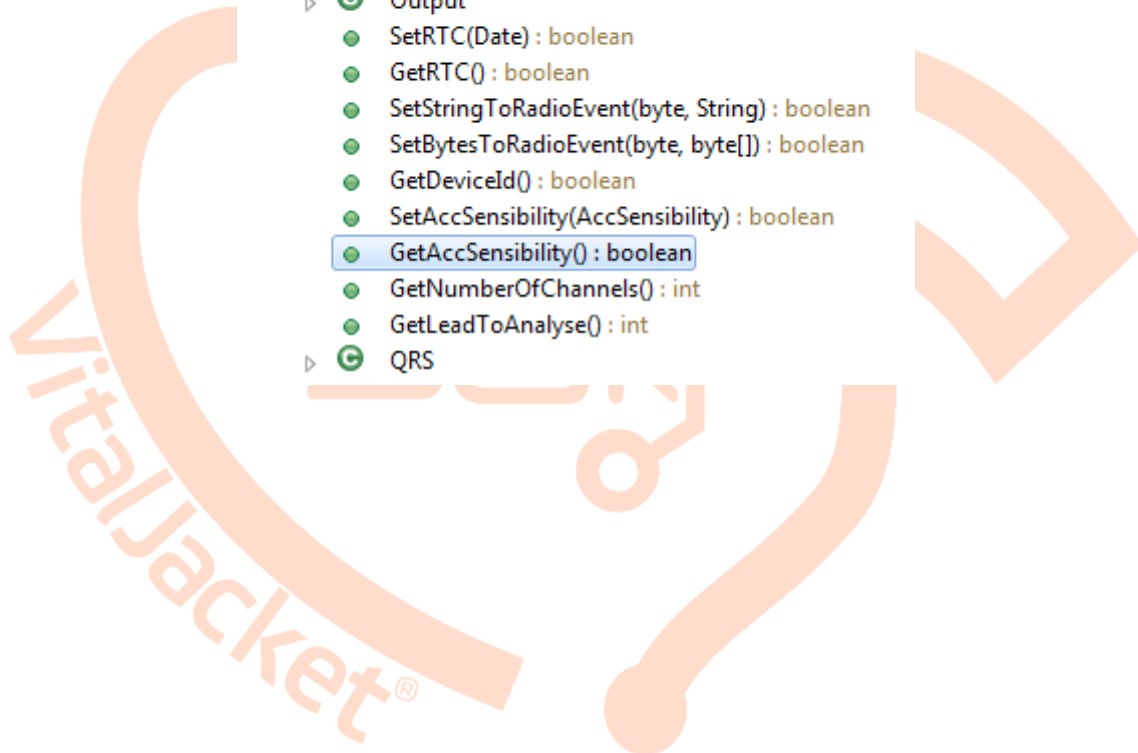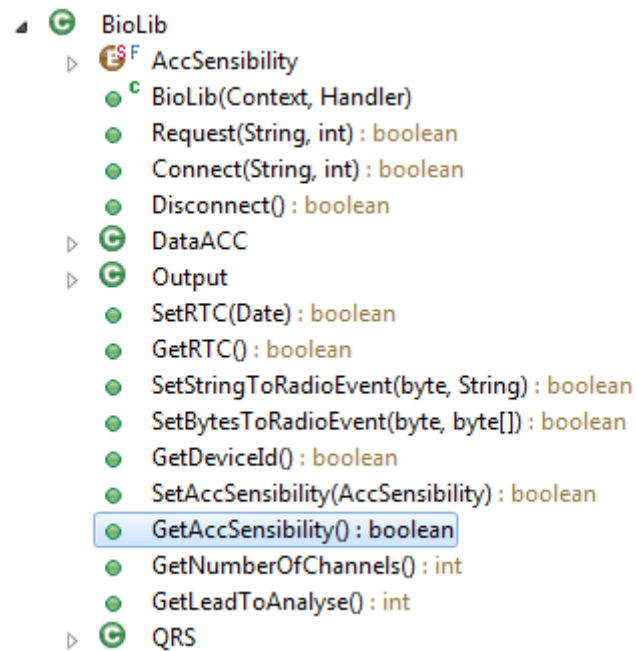

● **boolean Bio.Library.namespace.BioLib.Connect([String](#) address, int nQRS) throws [Exception](#)**

Connect to device (continuous mode). QRS detector is based on the algorithm of Pan and Tompkins [1] and was used MIT-BIH database for validate results [2].

Parameters:

address macaddress of device.

nQRS number of QRS to calculate BPM.

Returns:

true, if no errors occurred.

Throws:

[Exception](#)

```
BioLib lib = new BioLib(this, mHandler);

String address = "00:23:FE:00:0B:5E";

lib.Connect(address);
```

● **boolean Bio.Library.namespace.BioLib.Disconnect() throws [Exception](#)**

Disconnect from device (continuous mode).

Returns:

true, if no errors occurred.

Throws:

[Exception](#)

```
lib.Disconnect();
```

● **boolean Bio.Library.namespace.BioLib.Request([String](#) address, int time) throws [Exception](#)**

Connect to device (request mode).

Parameters:

address macaddress of device.

time duration of acquisition data.

Returns:

true, if no errors occurred.

Throws:

[Exception](#)

Note: Request for information to the device (heart-rate, push-button and battery level). Establishes a Bluetooth RFCOMM with the device, at the end of the preset time the connection is terminated.

```
lib.Request(address, 30);
```

● **boolean Bio.Library.namespace.BioLib.SetRTC([Date](#) date) throws [Exception](#)**

Set device RTC timecode.

Parameters:

date: date to update RTC.

Returns:

true, if no errors occurred.

Throws:

[Exception](#)

```
Date date = new Date();
lib.SetRTC(date);
```

● **boolean Bio.Library.namespace.BioLib.GetRTC() throws [Exception](#)**

Get RTC timecode from device.

Returns:

true, if no errors occurred.

Throws:

[Exception](#)

● **int Bio.Library.namespace.BioLib.GetNumberOfChannels()**

Get number of ECG channel(s).

Returns:

number of channel(s).

● **boolean Bio.Library.namespace.BioLib.GetDeviceId() throws [Exception](#)**

Get device Id from device.

Returns:

true, if no errors occurred.

Throws:

[Exception](#)

- **boolean Bio.Library.namespace.BioLib.GetFirmwareVersion()** throws [Exception](#)

Get firmware version from device VitalJacket.

Returns:

true, if no errors occurred.

Throws:

[Exception](#)

- **boolean Bio.Library.namespace.BioLib.SetBytesToRadioEvent(byte type, byte[] info)** throws [Exception](#)

Set radio event mark in ECG data.

Parameters:

type of radio event.

info data of radio event (10bytes maximum).

Returns:

true, if no errors occurred.

Throws:

[Exception](#)

```
// Sample of radio event
byte type = 1;
// Maximum 10 bytes to send device [Optional]
byte[] info = new byte[4];
info[0] = 0x31; // '1' ascii table
info[1] = 0x32; // '2' ascii table
info[2] = 0x33; // '3' ascii table
info[3] = 0x34; // '4' ascii table
lib.SetBytesToRadioEvent(type, info);
```

- **boolean Bio.Library.namespace.BioLib.SetStringToRadioEvent(byte type, string info)** throws [Exception](#)

Set radio event mark in ECG data.

Parameters:

type of radio event.

info data of radio event (10 chars maximum).

Returns:

true, if no errors occurred.

Throws:

[Exception](#)


```
// Sample of radio event
byte type = 2;
// Maximum 10 characters to send device [Optional]
String info = "5678";
lib.SetStringToRadioEvent(type, info);
```

● **boolean Bio.Library.namespace.BioLib.SetAccSensibility([AccSensibility](#) sensibility) throws [Exception](#)**

Set accelerometer sensibility (electronic device).

Parameters:

sensibility of accelerometer

Returns:

true, if no errors occurred.

Throws:

[Exception](#)

● **boolean Bio.Library.namespace.BioLib.GetAccSensibility() throws [Exception](#)**

Get accelerometer sensibility (electronic device).


Returns:

true, if no errors occurred.

Throws:

[Exception](#)

7. Print screens of Android application:

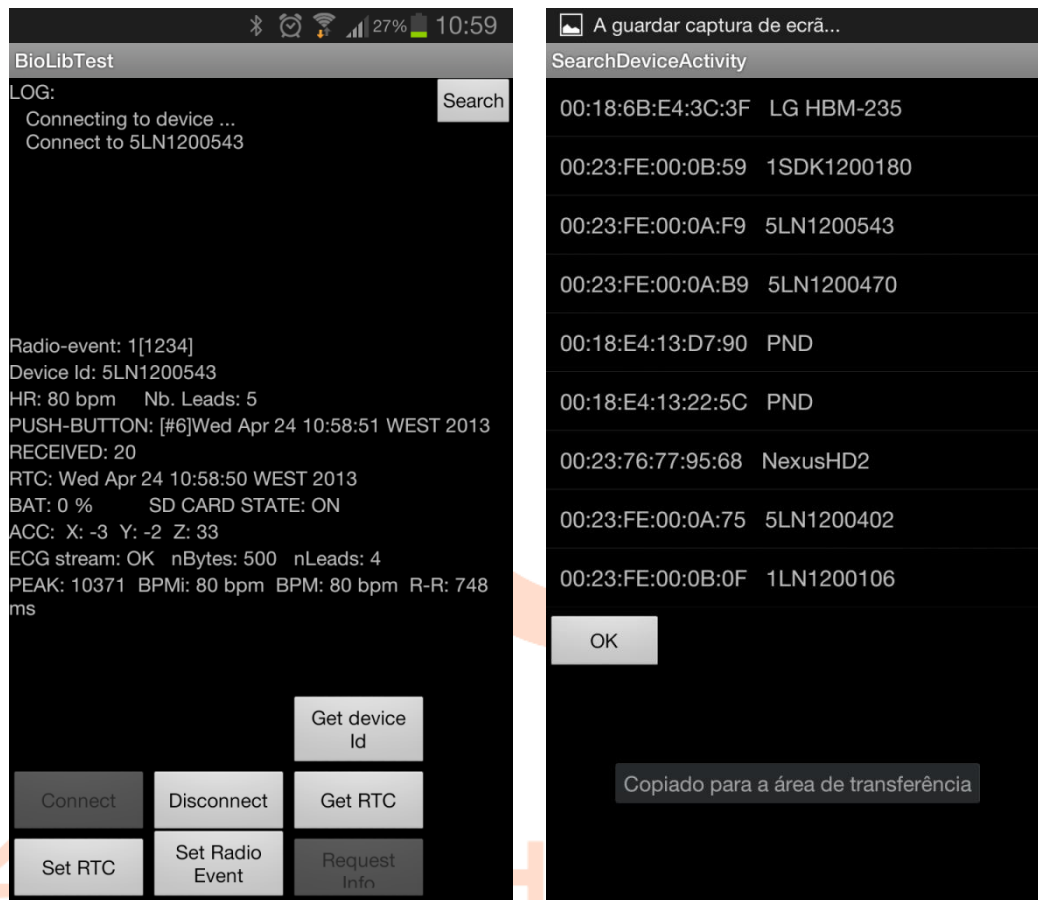

8. References

[1] Pan J and Tompkins WJ. A Real-Time QRS Detection Algorithm. IEEE Transactions on Biomedical Engineering 32(3):230-236, 1985

[2] MIT-BIH Arrhythmia Database: <http://www.physionet.org/physiobank/database/mitdb/>



9. Control versions

Version	Date	Change log
1.0.02	30-04-2013	<i>Get device ID</i>
		<i>Send radio event to device</i>
1.0.03	19-07-2013	<i>New method to send radio-event to device</i>
		<i>A new tool (InfoExporter.exe) for export data to Excel and Matlab</i>
1.0.04	01-05-2014	<i>Set / Get accelerometer sensibility</i>
1.0.07	18-03-2015	<i>Get firmware version from device VitalJacket.</i>

